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
EN PLEASE READ THIS OPERATING MANUAL CAREFULLY BEFORE STARTING THE ARTICLE. KEEP IT IN A SAFE PLACE FOR
FUTURE REFERENCE.

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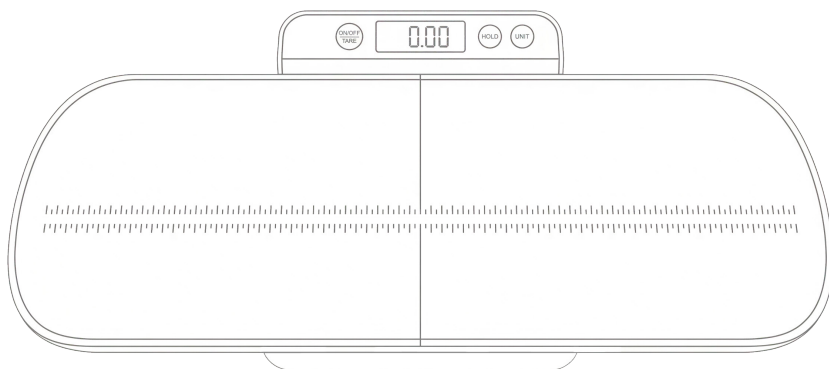EN THANK YOU FOR CHOOSING BELKO.

The BELKO Digital Baby Scale Plus enables precise and convenient measurements in everyday use. With features such as TARE, HOLD, and the integrated measuring ruler, you can reliably track your child's weight and length. Thanks to the removable tray, the scale can also be used as a regular body scale and is ideal for weighing pets.

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EN Do you have a question about the product? Encountered a problem? For technical support and service inquiries, you can confidently reach out to our trained staff.

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Please read these instructions carefully and follow all instructions, guidelines, and warnings included in this product manual in order to ensure that you install, use, and maintain the product properly at all times. These instructions MUST stay with this product. By using the product, you hereby confirm that you have read all instructions, guidelines, and warnings carefully and that you understand and agree to abide by the terms and conditions as set forth herein. You agree to use this product only for the intended purpose and application and in accordance with the instructions, guidelines, and warnings as set forth in this product manual as well as in accordance with all applicable laws and regulations. A failure to read and follow the instructions and warnings set forth herein may result in an injury to yourself and others, damage to your product or damage to other property in the vicinity. This product manual, including the instructions, guidelines, and warnings, and related documentation, may be subject to changes and updates. For up-to-date product information, please visit <https://belko-germany.de/manuals/>.

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EN SAFETY INSTRUCTIONS**WARNING!**

Please follow the safety instructions below to avoid injury or damage to the product. Read this manual carefully and completely before using the product for the first time. Keep this manual for future reference. If you pass the product on to another person, please also include this manual.

- Never leave your baby unattended on the scale.
- Always place the scale on a flat, stable, and solid surface.
- Do not use the scale on carpets or soft surfaces, as this may affect measurement accuracy.
- Do not use the product if it is damaged.
- Do not expose the scale to strong impacts, moisture, or water.
- Keep packaging materials and small parts away from children. Risk of suffocation!
- This scale is not a medical device and must not be used for medical diagnosis.
- Use the product only for its intended purpose.
- The product is intended for indoor household use only.
- Ensure that the tray is correctly assembled and locked in place before placing your baby on the scale.

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EN INTENDED USE

This digital baby scale is designed for measuring the body weight of babies, toddlers, or small pets. With the tray attached, the scale can be used for safely weighing babies. The measuring range in this mode is 0.1 kg to 40 kg.

After removing the tray, the scale can be used as a regular body scale.
The measuring range in this mode is 0.05 kg to 150 kg.

The product is intended for private household use indoors only.
For accurate measurements, the scale must always be placed on a flat, stable, and solid surface.

Improper use

The scale must not be used:

- for weights above 40 kg when the tray is attached
- for weights above 150 kg when the tray is removed
- as a step, seat, or toy
- on soft or uneven surfaces such as carpets or mattresses
- in humid environments or outdoors
- for medical or diagnostic purposes
- if the product shows visible damage


Using the product outside its intended purpose may lead to inaccurate measurements, product damage, or personal injury.

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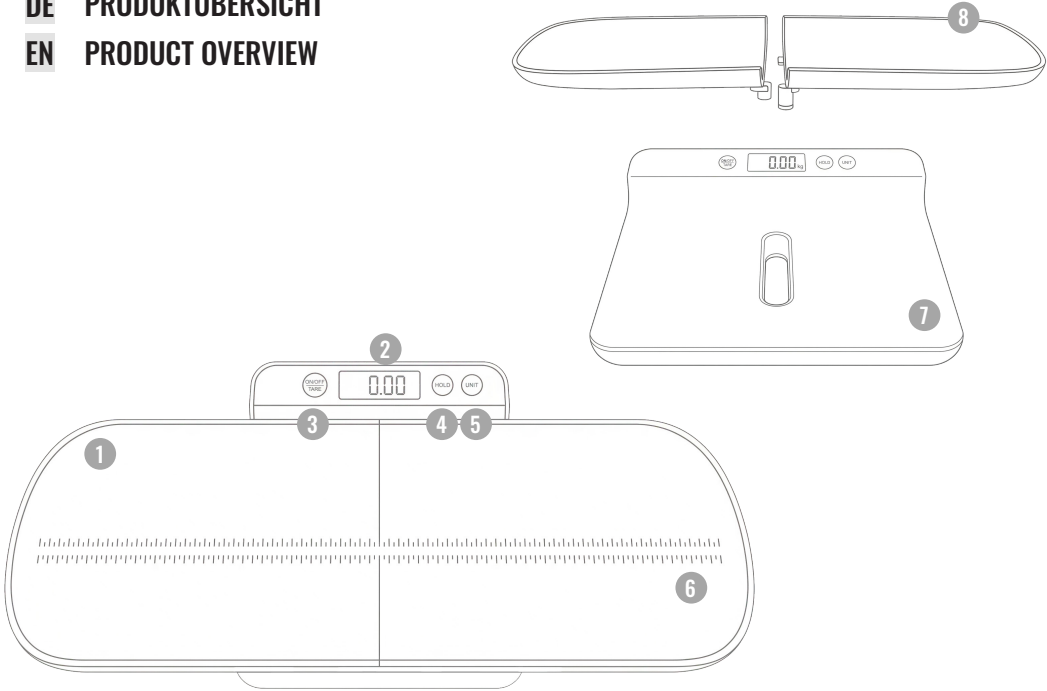
EN INSTALLING THE BATTERIES

1. Open the battery compartment on the underside of the scale.
2. Insert two AAA batteries.
3. Ensure the correct polarity (+ / -) when inserting the batteries.
4. Close the battery compartment securely.

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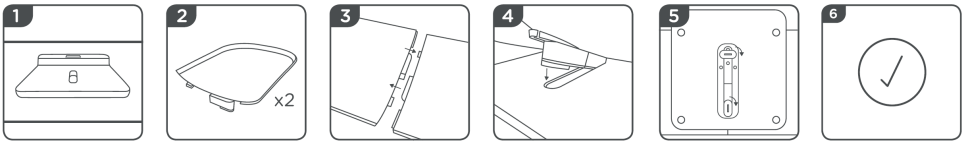
EN ASSEMBLY

1. Place the scale on a flat surface.
2. Remove the two baby tray sections from the packaging.
3. Push the two tray sections firmly together as shown.
4. Place the assembled baby tray onto the scale, aligning it with the mounting points
5. Turn the scale upside down and lock the tray in place by rotating the two locking discs 90°.
6. The assembly is now complete.

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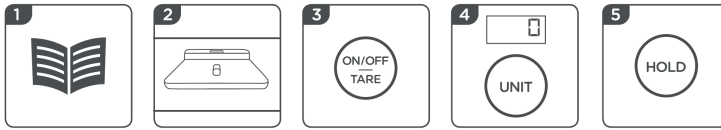EN **USING THE SCALE**

1. Place the scale on a flat and stable surface.
2. Press the ON/TARE button to turn on the device.
3. Wait until 0.00 appears on the display.
4. If using a blanket or towel, place it on the scale first and press TARE until 0.00 appears.
5. Carefully place the baby into the baby tray.
6. The weight will be shown on the display.
7. Use the UNIT button to change the weight unit (kg, oz, lb, lb:oz).
8. The HOLD function locks the weight on the display even if the baby moves during the measurement. Press the button again to cancel the function.

Note:

The scale can be used with the baby tray as a baby scale. Small pets can also be weighed in the tray. Remove the tray to use the scale as a regular personal scale.

The TARE function allows you to subtract the weight of items such as a blanket before weighing the baby.

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EN FUNCTIONS

TARE

Allows you to subtract the weight of blankets, towels, or clothing before weighing.

HOLD

Locks the measured weight on the display, even if the baby moves during the measurement.

UNIT

Changes the unit of measurement. Available units:
kg, oz, lb, lb:oz

Automatic Power-off

The scale will automatically switch off if it has not been used for 2 minutes.

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EN STORAGE

- When not in use, store the scale in a cool and dry place away from heat sources.
- If the scale will not be used for a long period, remove the batteries from the device.

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EN CLEANING AND MAINTENANCE

- Do not immerse the scale in water or other liquids, as this may damage the device.
- Clean the scale using a dry or slightly damp cloth only. Do not use harsh chemicals or cleaning agents.
- The baby tray can be washed with warm soapy water if necessary.
- Ensure all parts are completely dry before reassembly.

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EN TROUBLESHOOTING

Problem	Cause	Solution
Display does not turn on	Batteries are empty or inserted incorrectly	Replace the batteries and check the correct polarity (+/-)
“LO” appears on the display	Battery low	Replace with new AAA batteries
“ERR” appears on the display	Maximum weight exceeded	Remove the weight immediately and do not exceed the maximum capacity
Incorrect weight reading	Scale is placed on carpet or soft surface	Place the scale on a hard, flat surface such as wood or tile flooring
Display does not show “0.000” before weighing	Scale was not reset before weighing	Press the ON/TARE button to reset the scale
“88888” appears on the display	Device is restarting	Reset the battery or wait until the display jumps to „0“
- (before the digits)	Imbalance or negative value	Place the scale on a hard, flat surface and press TARE
“Un5L” appears on the display	Imbalance	Place the scale on a hard and flat horizontal surface, and with no other objects touching the scale.

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EN TECHNICAL SPECIFICATIONS

Specification	Digital Baby Scale
Maximum capacity (with tray) Maximum capacity (without tray)	40 kg 150 kg
Measurement range (with tray) Measurement range (without tray)	0.1 kg - 40 kg 0.05 kg - 150 kg
Measurement accuracy	0.05kg - 4.99kg: ±50g 5.00kg - 19.99kg: ±100g 20.0kg - 49.0kg: ±200g 50.0kg - 99.9kg: ±300g 100.0kg - 150.0kg: ±400g

EN TECHNICAL SPECIFICATIONS


Specification	Digital Baby Scale
Power supply	2 × AAA batteries
Product weight	2 kg
Product Dimensions	73.2 x 34.5 x 6.5 cm
Single tray section size Assembled baby tray size	38 x 27 cm 73,6 x 27 cm
Height measuring scale	0 - 69 cm
Package Dimensions	416 × 337 × 83 mm

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
EN BATTERY REPLACEMENT

- Batteries should only be installed or replaced by an adult.
- Use only the specified battery type (2 × AAA, 1.5 V) and install them correctly.
- Do not mix different battery types or new and used batteries.
- Keep batteries out of reach of children.
- Dispose of used batteries in accordance with local regulations.
- Do not dispose of batteries in fire.

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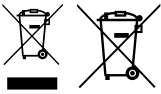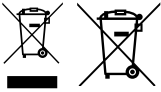EN DISPOSAL

Electrical devices must not be disposed of with household waste. This product must be disposed of in accordance with the European WEEE Directive (Waste Electrical and Electronic Equipment). Please dispose of this product at an appropriate collection point for the recycling of electrical and electronic equipment.

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EN WARRANTY

If the device becomes defective, a 1-year warranty from the date of purchase is granted for this item, except in the following cases:

- Willful damage
- Damage caused by force majeure, such as flooding, fire, etc.
- Damage due to improper use or disregard of these instructions.
- Damage or malfunction caused by disassembling or modifying the device.

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EN EU DECLARATION OF CONFORMITY

The manufacturer / entity placing the product on the market
Panorama24 GmbH
Chromstr. 4
30916 Isernhagen, Germany
Tel. (+49) (0) 511 - 21241670
hereby declares that the following product

Product designation: Belko® Digital Baby Scale Plus
Model designation:

P24BSCALEPLUS (B2H, Belko® Digital Baby Scale Plus)

Year of manufacture: 2026

complies with all the relevant provisions of the applied legislation (listed below) — including any amendments valid at the time of this declaration. The manufacturer assumes sole responsibility for issuing this Declaration of Conformity. This declaration applies only to the machine in the condition in which it was placed on the market; any parts added or modifications made by the end user at a later date are not taken into account.

The following legislation was applied:

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EN APP FUNCTION

The “BabyFit” app is a mobile application designed to track and analyze your baby’s development. It records measurement data such as weight, height, head circumference (HC) and BMI development. All measurements are clearly displayed, allowing you to easily monitor your baby’s growth. The app supports up to 8 user profiles.

Scan the QR code matching your operating system or search for the „BabyFit“ app in the App Store for Apple devices or the Google Play Store for Android devices.

Step 1: Install the App

- Scan the QR code or download the “BabyFit” app from the App Store or Google Play Store.
- Open the app and create a user account.

Step 2: Connect the scale

- Enable Bluetooth on your smartphone.
- Turn on the scale.
- Open the app and add a new device.
- Select the baby scale from the device list.

Step 3: Create a baby profile

Select the desired weight unit and confirm by tapping “OK”.
Then Enter your baby’s information, such as:

- name
- gender
- date of birth, then tap “Confirm”.

Step 4: Perform a measurement

- Place the baby on the scale.
- The weight will automatically sync with the app.
- The app allows you to view measurement history and growth charts.

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Copyright 2026 Panorama24 GmbH
Belko® ist eine Marke der Panorama24 GmbH